
YOX

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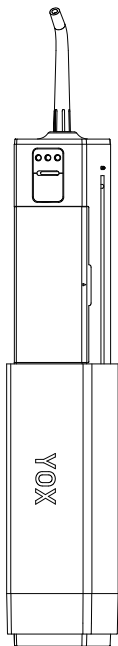
Water Flosser
Model: YJW-600

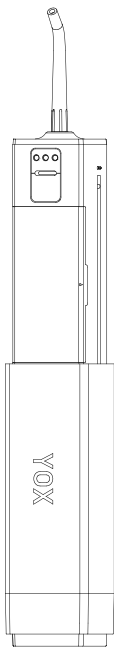Before Use

Thank you for purchasing the YOX Water Flosser. Please read this user manual carefully and use the product correctly and safely.

Safety Instructions

Warning

Stop using the product immediately if there's any abnormality.

Keep out of reach of infants and young children. Do not allow them to touch it.

Note

Prohibited Behaviors (Things you should not do):

Do not drop, step on or hit it.

Do not leave it in a car.

Do not expose it to direct sunlight.

Do not place it near heat sources such as cooking utensils or heating devices.

Do not overload it with high temperature, shock, or overcharging.

Do not modify or disassemble it.

Do not remove the charging port cover.

Do not put water into the charging port.

Do not use it if damaged.

Do not use it with alcohol or corrosive liquids (not covered by warranty).

Do not use it with warm water (up to 40°C).

Do not allow people who cannot operate it themselves to use it.

Do not allow people who cannot communicate or have physical disabilities to use it.

Do not use this product on persons with reduced mobility.

Do not use it for purposes other than oral cleaning.

Do not share the tip with others.

Do not charge during a short circuit.

Do not charge it with anything other than the included special USB cable.

Do not use the accessories for any other purpose than this product.

Note

Compliance items (please follow):

Always keep the product clean.

After use, drain the water from the reservoir and let it dry naturally.

Clean the tip and reservoir frequently, and let them dry completely.

If the user is a child, elderly person, or someone who has never used the product before, the use should be agreed, supervised and instructed by a parent or guardian who is familiar with the product.

If you have symptoms of periodontal disease or other oral problems, consult a doctor before using the product.

If you use mouthwash, use one that does not contain alcohol.


Turn off the power after use.

Do not fill the reservoir with water while charging.

Drain the water from the reservoir when not in use.

Dispose of the product as specified by the local authorities.

Names of Each Part

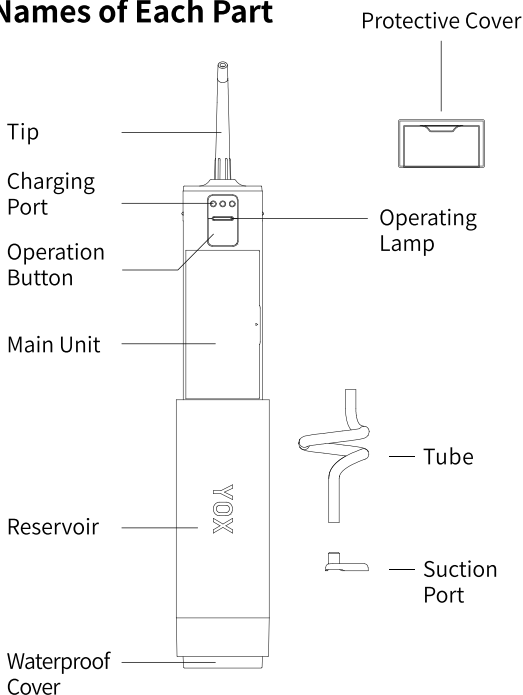The suction port is magnetic and adheres to the bottom of the reservoir. If it comes off and floats, reattach it to the bottom of the reservoir before use.

If the tube comes off the suction port, reinsert it before use.

Product Specifications

Product Name: YOX Water Flosser

Model	YJW-600
Dimensions	36 × 36 × 123mm
Weight	128g
Reservoir Capacity	72ml
Charging Time	Approximately 120 minutes※1
Continuous Usage Time	Approximately 60 minutes※1
Charging Port	Exclusive USB cable (Magnetic type)
Modes	Soft cleaning, strong cleaning
Mouthwash	Available for use※2
Specifications	5V/1A
Battery	500mAh
Waterproof Level	IPX7

※1

Charging time and continuous use time are estimates and are not guaranteed. They may vary depending on usage and environmental conditions.

※2

This applies to non-alcoholic mouthwash only.

Warranty Certificate

This is required when receiving warranty service in your region.
Please keep it safe.

Product Name	
Manufacturing Number	
Customer Information	Address/Name/Phone Number The personal information of the customer written on the warranty certificate will be used only for repairing or replacing the product. It will not be used for any other purposes.
Date of Purchase	Year/Month/Day
Warranty Period	Main Unit: 1 year from the date of purchase
Retailer	

YOX


Hm Brand Operation Management Limited
Email: info@yox.vip

Instructions

[1] Pull out the main unit

Open the protective cover and pull out the main unit from the reservoir.

Please be careful not to let the main unit fall out when pulling it out.


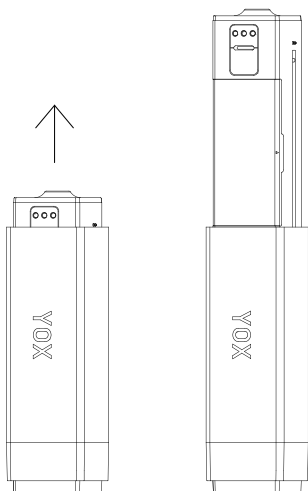
- ! This product is USB rechargeable. Please fully charge it before use. The operation lamp will blink during charging and will stay on constantly when fully charged. Please charge it with an empty reservoir.

[2] Attach the tip

Attach the tip to the top of the main unit.

While checking the position of the tip groove and the shape of the insertion port, insert the tip and rotate it 90 degrees to attach it in the correct position. If the front of the operation button and the tip are facing the same direction, it is properly attached.

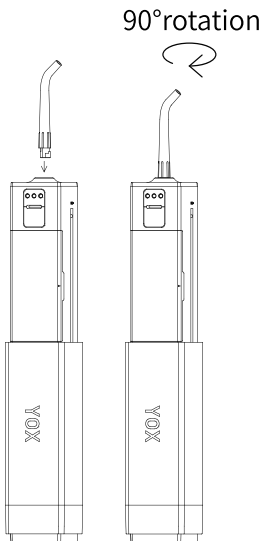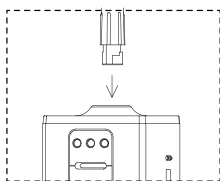
After inserting the tip into the main unit, rotate it about 90 degrees to lock it in place.

[3] Add water

Add water from the water supply port.

Turn the main unit upside down, slide the waterproof cover, and add water from the water supply port.

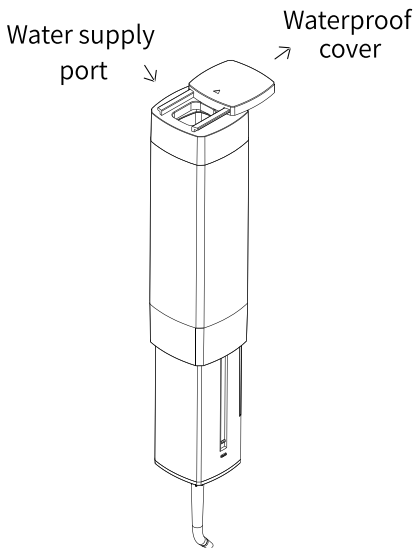[4] Oral cleaning

With your mouth closed, press and hold the operation button to start.

Align the position where you want to floss with the position of the tip, close your mouth, and then press and hold the operation button. It will start in Strong Cleaning mode.

Pressing the operation button once while in Strong Cleaning mode will switch it to Soft Cleaning mode. Pressing it once more will stop the operation.

If the water in the reservoir runs out, return to step 3 and add more water.

It will automatically stop after about 40 seconds of continuous operation, but you can use it again by operating it again.

Cautions

Pressing the operation button for a long time will cause water to spray out vigorously. After aligning the tip, be sure to close your mouth before operating the button.

When cleaning your mouth, water may spill out of your lips if you are not used to it. Please prepare a towel or other things before operating.

Depending on the operating method, there is a possibility that water may splash. Please refrain from operating near equipment that is sensitive to water.

The tip is a consumable.

If the tip is hard to assemble because the jet power is weaker, water leaks, or damaged, please replace a new tip.

Replace a tip every 3 months.

[5] Storage


Remove any remaining water and store

To store the product, remove any remaining water in the reservoir by sliding the water supply cover. Then, use a dry cloth or paper towel to wipe any remaining moisture from the product.

Remove the tip and store it in the tip storage compartment.

With natural ventilation, leave the water supply cover open.

Store the main unit in the clear case and attach the protective cover.

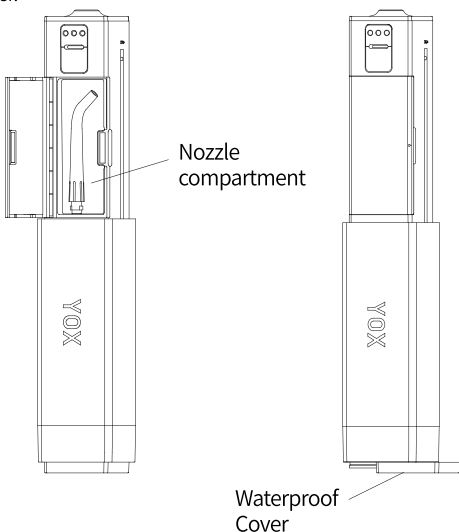
Charging:

Connect the charging cable and the charging port through a magnet.

Connect them in the right direction. When close to the wrong direction, the magnet will rebound. It will make the charging cable hard to reach the charging port. Please change the direction for connecting.

When connecting, the charging cable and the charge port may be deviated, but it can be used directly.

When charging, please use the original USB cable(magnetic)

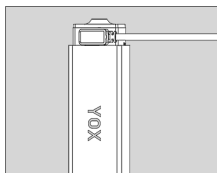
▲ Magnet connecting direction

Even if during the connecting, the charging cables and the charging port may be a little deviated.


※Cleaning the reservoir

Fill the reservoir with water.

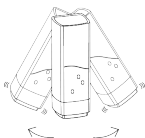Shake the reservoir from side to side.

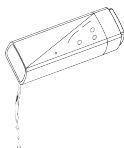
Discard the water.

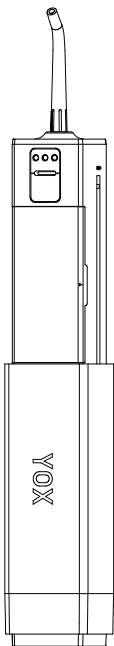Allow to dry naturally. During drying, partially open the water supply cover to improve ventilation.

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